

Introducing

# ErbaQik<sup>®</sup>

One-step Rapid Tests  
for **Critical Infections**

**ErbaQik<sup>®</sup> HIV 1+2**

**ErbaQik<sup>®</sup> HBsAg**

**ErbaQik<sup>®</sup> HCV**

**ErbaQik<sup>®</sup> Syphilis**



**ErbaQik<sup>®</sup> HIV 1+2**



**ErbaQik<sup>®</sup> HBsAg**



**ErbaQik<sup>®</sup> HCV**



**ErbaQik<sup>®</sup> Syphilis**

## For Quick and Accurate Testing:

- Lateral flow Immunochromatographic (ICT) assay
- Sample type: Whole blood\*/ plasma/ serum
- Results in 20 mins
- Excellent sensitivity and specificity
- 24 months shelf life
- Storage temperature: 2-30°C


Evaluated by



National Institute of Biologicals (NIB)


\*Except for HBsAg

Transasia's ErbaQik rapid tests for critical infections are useful tools for screening of HIV<sub>1+2</sub>, HBsAg, HCV and Syphilis.

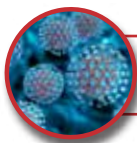

### ErbaQik<sup>®</sup> HIV<sub>1+2</sub> Test Card

- Detection of total antibodies (IgG, IgM & IgA)
- Differential detection of HIV-1&2 antibodies
- Test region 1 coated with recombinant antigen gp120 and gp41 for HIV-1 and test region 2 coated with gp36 for HIV-2
- Sensitivity: 100%, Specificity: 100%

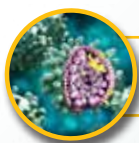


### ErbaQik<sup>®</sup> HBsAg Test Card

- Detection of HBsAg of all 11 subtypes and known mutant strains
- Test region coated with monoclonal anti-HBsAg antibodies
- Sensitivity: 100%, Specificity: 100%



### ErbaQik<sup>®</sup> HCV Test Card

- Detection of anti-HCV antibodies
- Test region coated with mixture of chimeric recombinant HCV antigens (Core, NS3, NS4 and NS5)
- Sensitivity: 100%, Specificity: 100%

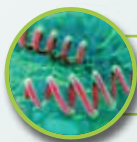

### ErbaQik<sup>®</sup> Syphilis Test Card

- Detection of antibodies against Treponema pallidum
- Test region coated with recombinant Treponemal antigen (Tpn15, Tpn47 and Tpn17)
- Sensitivity: 100%, Specificity: 100%

## Ordering information

Product Code	Product Name	Pack Size
132821	ErbaQik HIV 1+2 Test Card	25T
133156	ErbaQik HBsAg Test Card	25T
132825	ErbaQik HCV Test Card	25T
133160	ErbaQik Syphilis Test Card	25T